

**PRODUCT BROCHURE**

---

**KERN  
SET-37**

Date: 18.03.2025

**ENGLISH**

Software EasyTouch

SET-37

**Export  
Manager**

EasyTouch Export Manager – export results data automatically and reliably from EasyTouch – e.g. to import into your ERP system

---

**CONTACT****KERN & SOHN GmbH**

Ziegelei 1

72336 Balingen






Germany

Phone : +49 7433 9933-0

Fax : +49 7433 9933-149

Email : [info@kern-sohn.com](mailto:info@kern-sohn.com)Web : [www.kern-sohn.com](http://www.kern-sohn.com)

## PRODUCT SHOWCASE

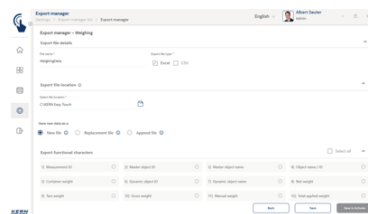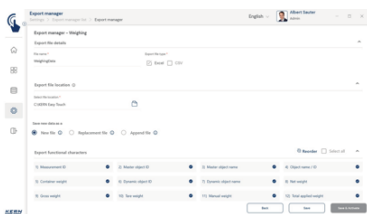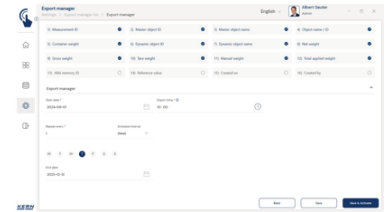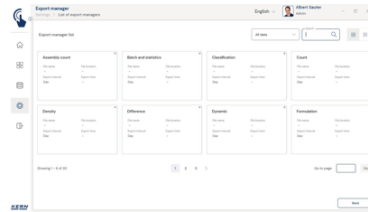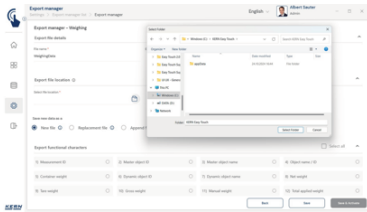

ENGLISH

## DESCRIPTION

- Note: Prerequisite for this set is the basic program SET-01 Base
- Data export: Export Manager offers automatic data export at a time interval which can be selected individually. When doing this, different file formats such as CDV or XLSX are supported, in order to enable seamless further processing of data, for example in ERP and PPS systems.
- Data storage location: The storage location for exported files can be freely defined
- Data fields: Depending on the EasyTouch function selected, the data fields to be exported can be freely defined.
- Export intervals: Times and intervals for exports can be set to suit, e.g. once or several times a day, weekly, monthly and annually at defined times
- Export history: Using Export log history you can access all previously exported files
- SaveServer: Use this function to export data from the central SaveServer
- Export file variants: You can choose between different export options:
  - a. Export to a new file
  - b. Export to an existing file

## SPECIFICATIONS

### Article codes, product groups & Status information

Model article-code	SET-37
Model-Serie	ETS
Assortment	KERN
Product Type	Software
Product group	EasyTouch
Product usage type	Main item;Accessories/Services
EAN code (Model)	404576 1476523
Customs tariff number	85234910

### Technical data - Packing & Shipping

Delivery time	1 d
---------------	-----

## FUNCTIONS

Standard



Option

