
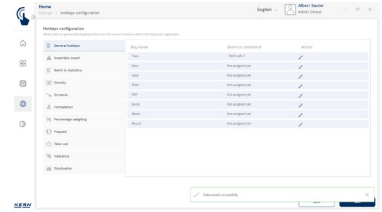
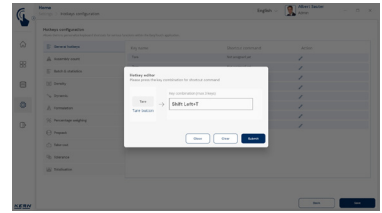


Software EasyTouch SET-411 Hotkeys

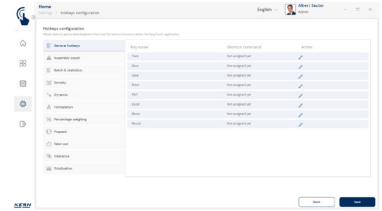

EasyTouch hotkeys - shortcut keys and external switches

Features

- Note: Prerequisite for this set is the basic program SET-01 Base
- Shortcut keys: Hotkeys allow you to assign key combinations to all buttons in EasyTouch. This simplifies access to frequently used functions considerably, as it means the laborious navigation with the PC mouse is not required. Instead, the actions required can be triggered quickly and directly using the keyboard
- Optional item: Hotkeys is an optional extension of EasyTouch, which supports all active, installed functions. Please note that the price for the additional function does not include any hardware such as switches, buttons or foot pedals
- Configuring shortcut keys: Key combinations can be defined by the system administrator or an authorised person. This means that a clear, user-specific key combination can be defined for any required function
- Versatile hotkey selection: In EasyTouch, you can define general key commands such as "Save" or "Tare" as hotkeys, as well as specific action keys for certain functions. These can then be used in all EasyTouch functions
- Activation using key combination: For example the TARE function can be activated simply and quickly using a key combination
- Integration of external switches or keys: EasyTouch supports the assignment of key combinations to signals from external HID devices (Human Interface Devices) via USB. This allows the software to be operated via external switches such as emergency stop buttons or foot switches, which enables hands-free operation and increased efficiency. About the process: The external signal is assigned to a key combination in the control software of the standard USB device. This is received by the PC. The signal then activates the assigned EasyTouch function

