
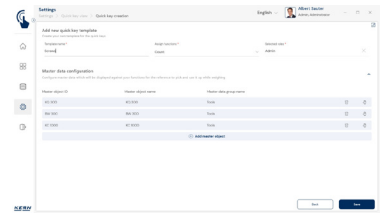


Software EasyTouch

SET-311

Quick Key

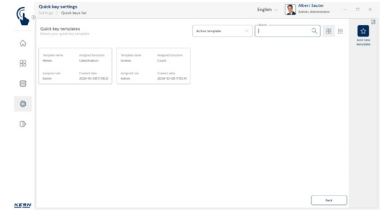

EasyTouch Quick Key – select master data objects quickly and user-related

Features

- Prerequisite for this set is the basic program SET-01 Base
- Quick Keys can be created in numerous sets, e. g.
 - SET-01 Base
 - SET-05 Dynamic
 - SET-06 Classification
 - SET-07 Variable
 - SET-08 Take Out
 - SET-31 Tolerance
 - SET-32 Count
 - SET-33 Target Count
- Quick selection keys for master data: Through ET Quick Keys, specific master data objects can be permanently displayed for quick selection. The Quick Keys are displayed on the screen immediately

- Saving: Whether selected in the conventional way or by using Quick Keys, every result of a weighing or piece count will be stored in relation to the selected master data object
- Selecting master data objects: The specific master data objects, which are displayed as Quick Keys, can either be selected individually or as a defined master data group, which contains the master data as Quick Keys
- User-related: The Quick Keys can be created as a scheme which relates to particular user roles. This means that the users of a particular role are shown the relevant master data objects for quick selection. This allows, for example, different departments in the same company to work with different master data without overlapping

License information

- User: You can store as many users as you need in one license
- Balances: You can store and operate as many balances as you need in one license
- Communication between balance/terminal device: The balance(s) can communicate with the PC, laptop or tablet by serial connection, USB, Bluetooth, Ethernet or WiFi

STANDARD



OPTION

