KERN & Sohn GmbH
Ziegelei 1
72336 Balingen-Frommern
Germany
+0049-[0]7433-9933-0
+0049-[0]7433-9933-149
info@kern-sohn.com

# Installation Instructions Stands

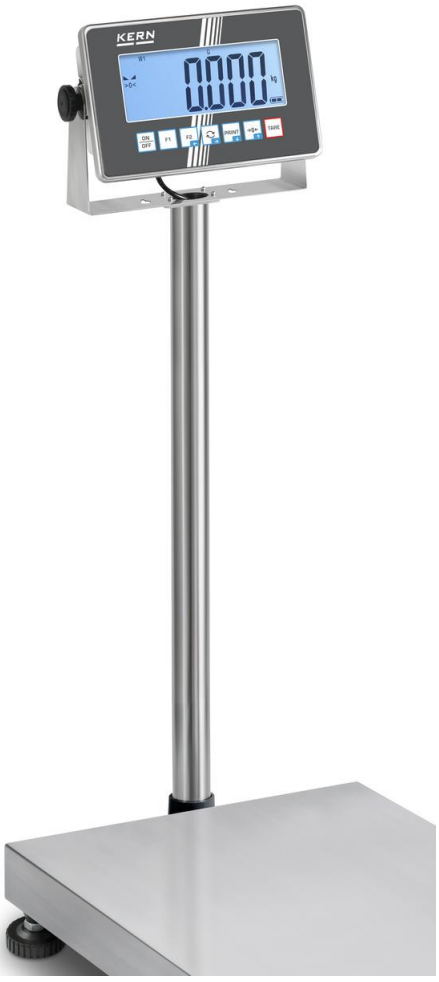
## KERN SXC, IXC

Type

TSXC-A01-A
TIXC-A01-A
TIXC-A02-A

Version 2.0
2026-06
en





en

# Stands

Version 2.0 2026-06

## Installation Instructions

Contents


**1 Technical data ..... 3**

**1.1 Dimensions ..... 3**

**2 Scope of delivery ..... 6**

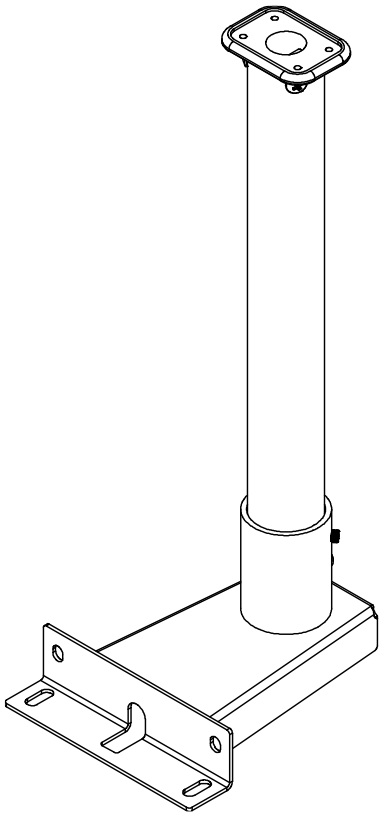
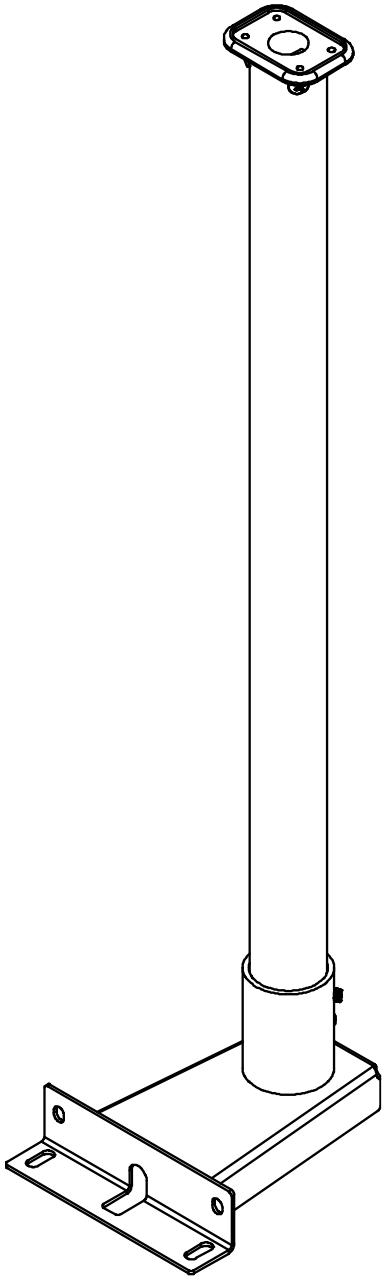
**3 Installation ..... 7**

# 1 Technical data

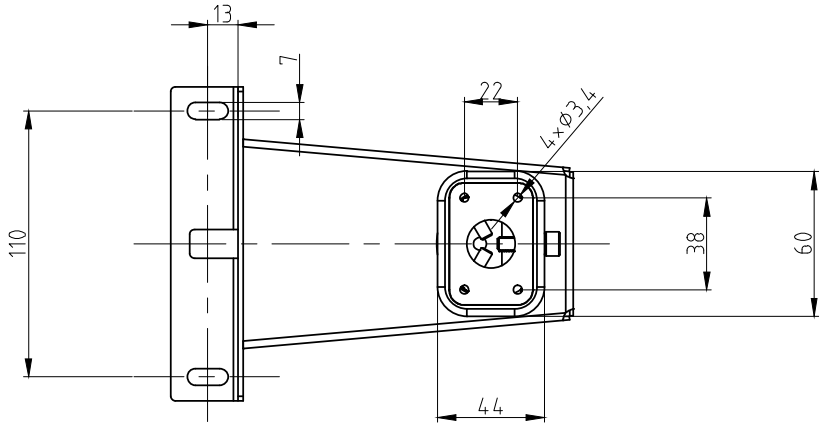
| KERN | TSXC-A01-A | TIXC-A01-A | TIXC-A02-A |
|---|---|---|---|
| Stand height (mm) | 330 | 330 | 600 |

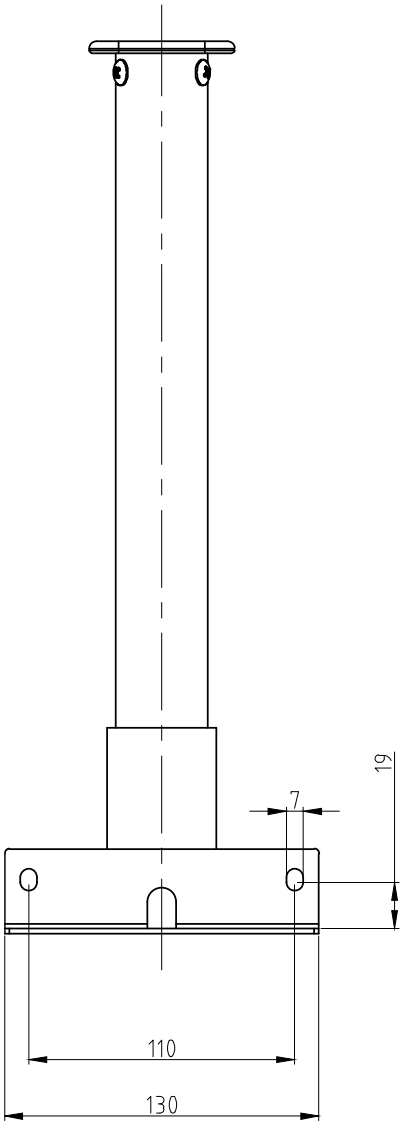
## 1.1 Dimensions

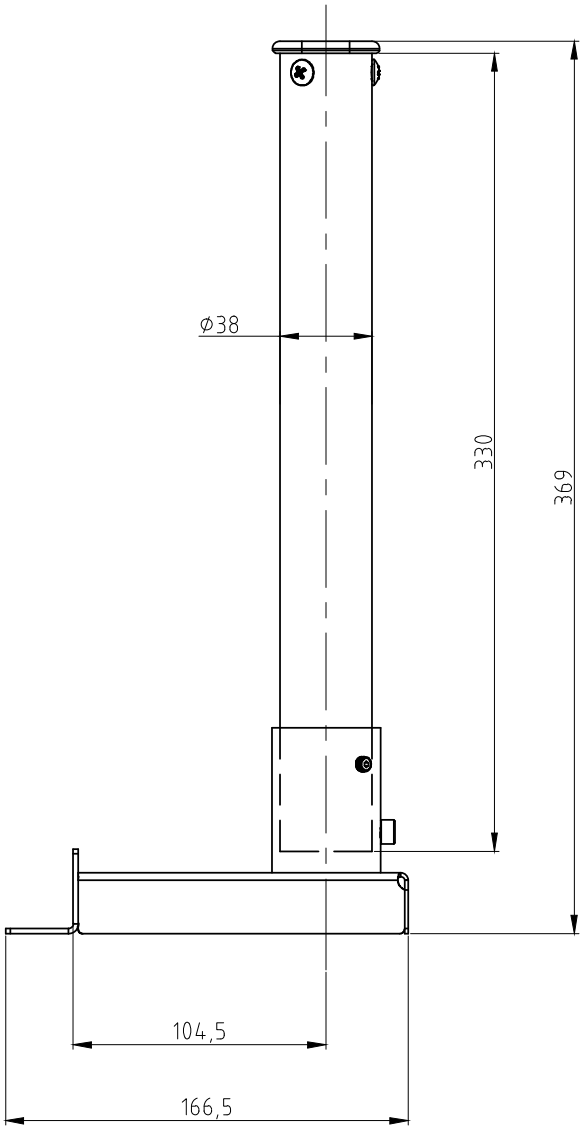
- Dimensions in mm
- Drawings on the following page

| Models with 330 mm | Models with 600 mm |
|---|---|
|  |  |

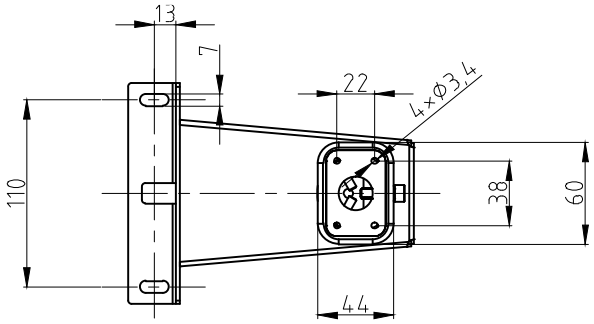
**Models with 330 mm:**

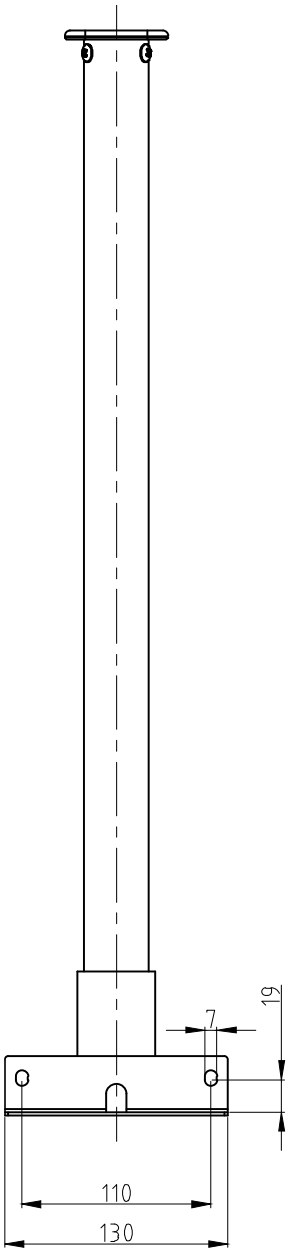


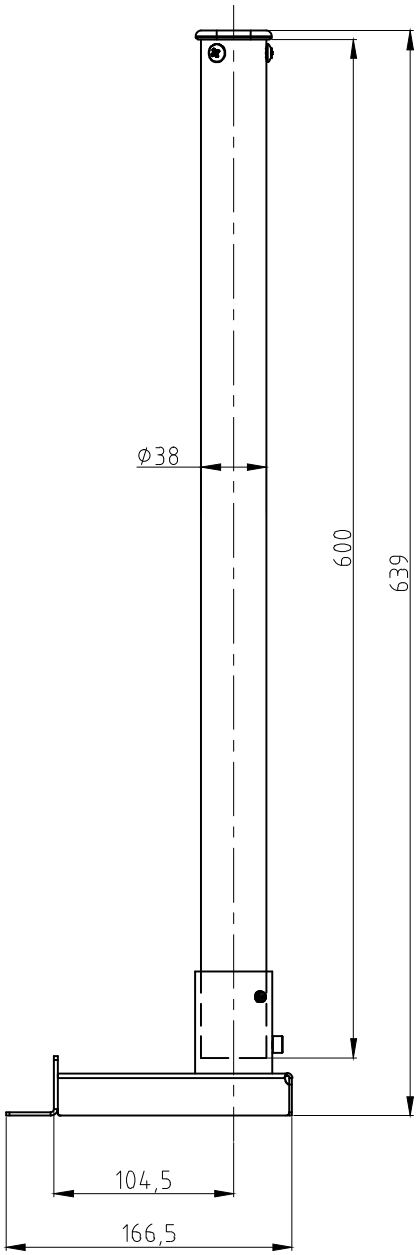



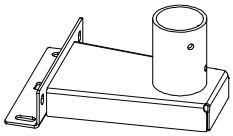
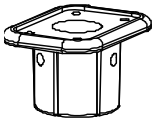
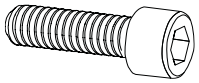
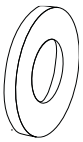
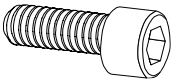
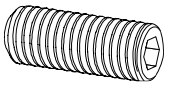
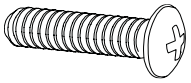
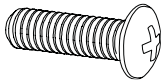
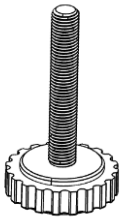
**Models with 600 mm:**



## 2 Scope of delivery

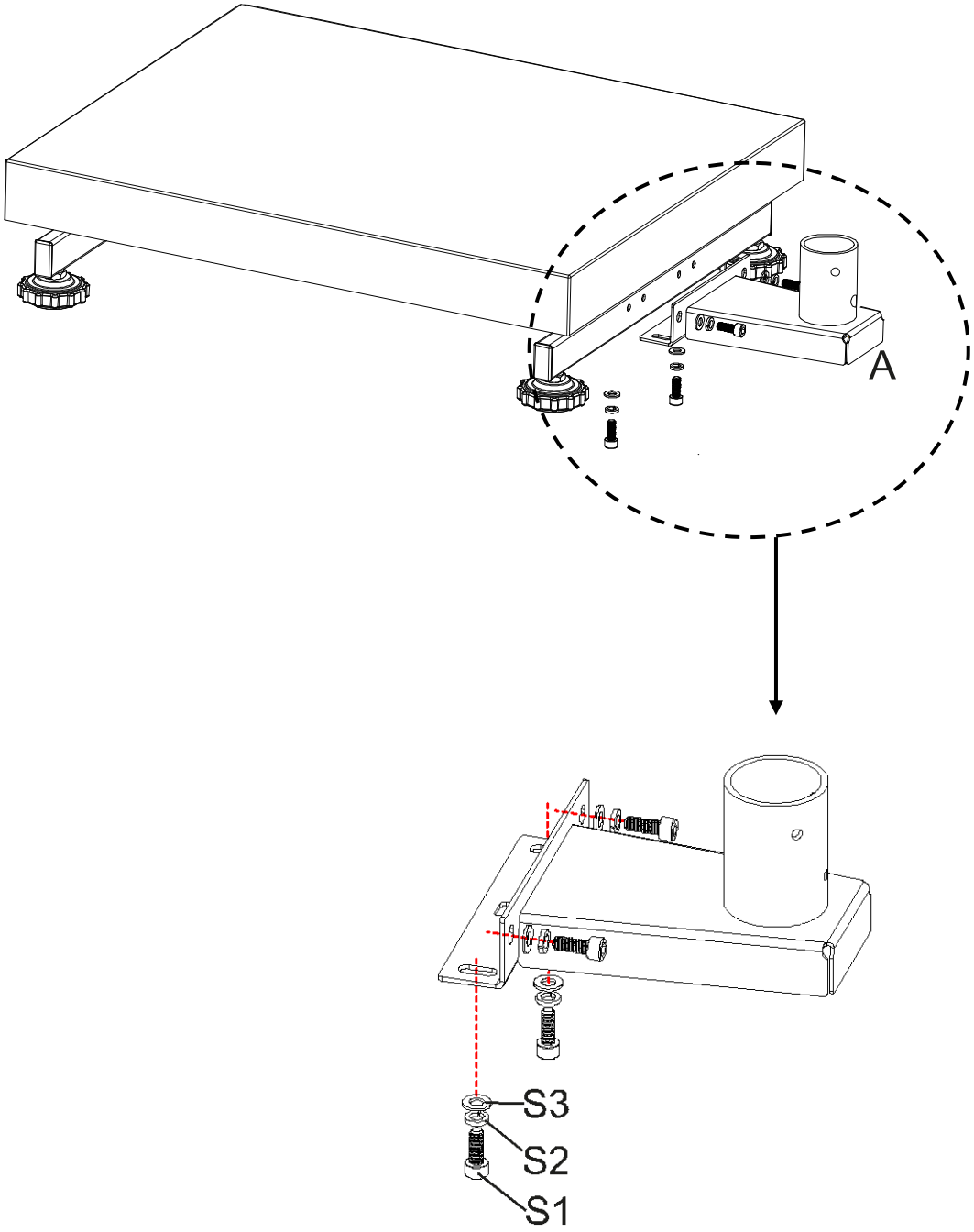
| | | |
|---|---|---|
| <br>**A**<br>Stand foot<br>1 x | <br>**B**<br>Stand tube<br>1 x | <br>**C**<br>Support for display unit<br>1 x |
| <br>**S1**<br>M6 x 16<br>4 x | <br>**S2**<br>Spring washer<br>4 x | <br>**S3**<br>Washer<br>4 x |
| <br>**S4**<br>M6 x 12<br>1 x | <br>**S5**<br>M6 x 10<br>2 x | <br>**S6**<br>M5 x 16<br>3 x |
| <br>**S7**<br>M4 x 10<br>4 x | <br>**S8**<br>Adjustable foot<br>1 x | <br>**S9**<br>M6 nut<br>1 x |

# 3 Installation

**Assembling the stand foot to the platform:**

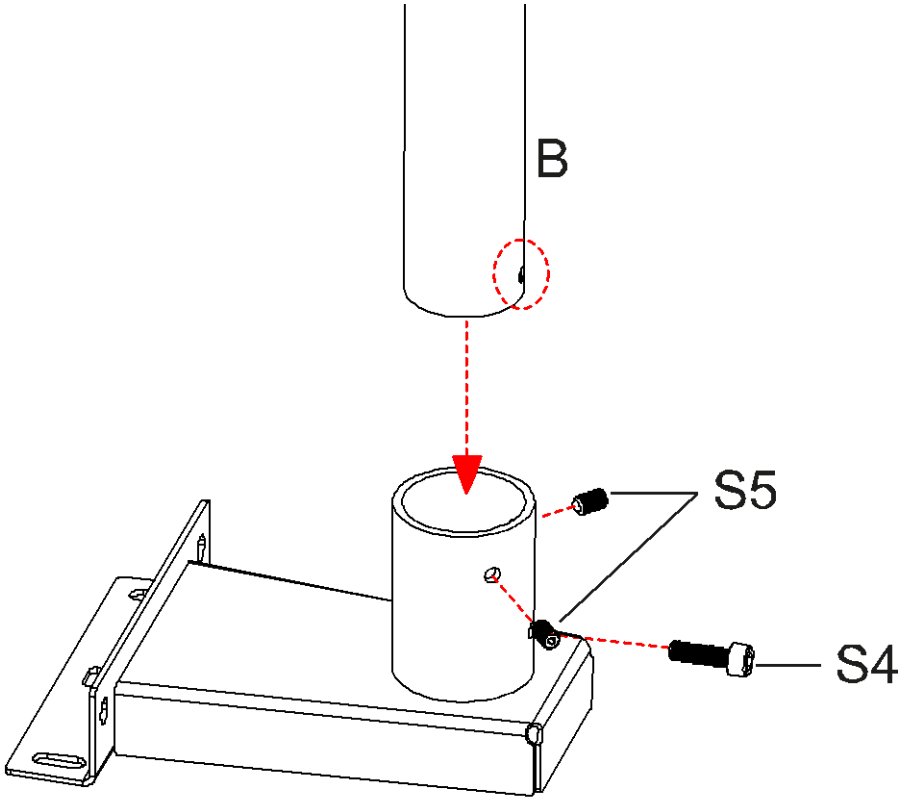
1. Place the stand foot **(A)** centrically over the perforations on the platform side.
2. First plug the spring washer **(S2)** and then the washer **(S3)** on the screw **(S1)** (4 x).
3. Screw down the stand foot.

   **Make sure that no cables are jammed or damaged.**

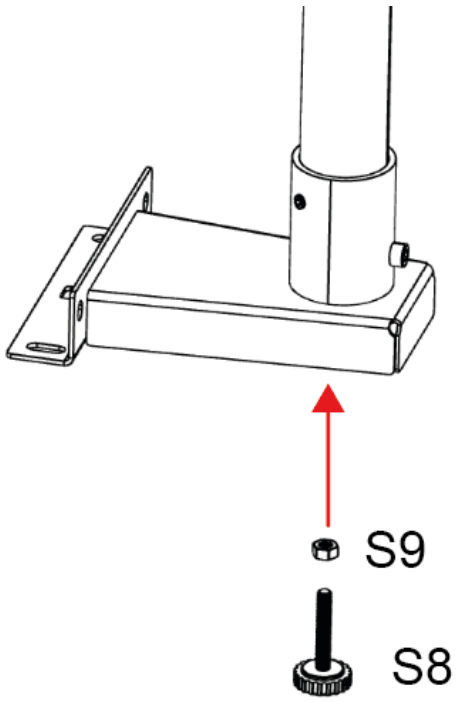


**Assembling the stand tube to the stand foot:**

4. Unscrew the screws **S4** and **S5** from the stand tube.
5. Push the stand tube **(B)** into the stand foot. Ensure that the perforation of the stand tube matches with the perforation **S4**.
6. Screw down again the screws **S5** and fasten the stand tube using screw **S4**.



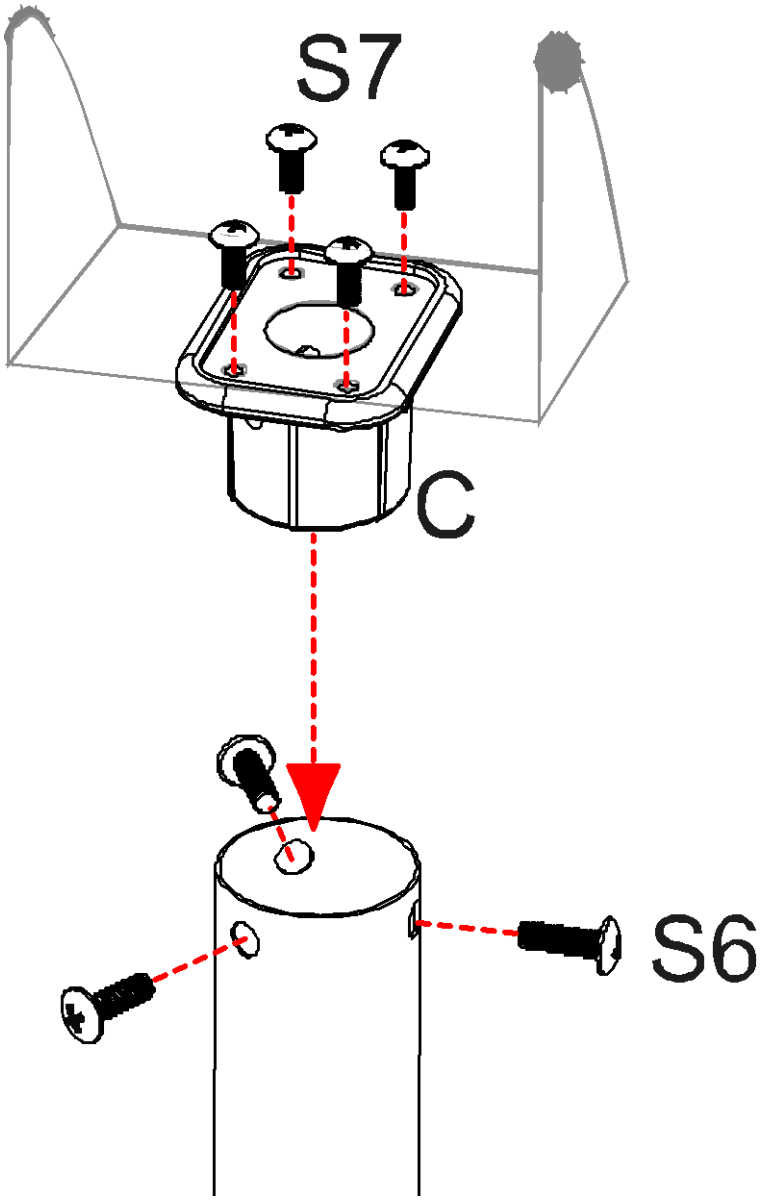

7. Screw nut **S9** onto the adjustable foot **S8**.
8. Screw the adjustable foot, with the nut, into the tripod tube from below.
9. Unscrew the adjustable foot until it touches the floor and secure it with nut **S9**.




**Assembling the display unit to the stand tube:**

**10.** Push the support of the display unit **(C)** into the stand tube and fasten it using screws **S6**.

**11.** Fasten the display unit to the display unit's support using screws **S7**.


S7
C
S6