

PRODUCT BROCHURE

**KERN
SET-40**

Date: 18.03.2025

ENGLISH

Software EasyTouch

SET-40

Protect

EasyTouch Protect – protect sensitive company data safely and efficiently

CONTACT

KERN & SOHN GmbH

Ziegelei 1

72336 Balingen

Germany



Phone : +49 7433 9933-0

Fax : +49 7433 9933-149

Email : info@kern-sohn.com

Web : www.kern-sohn.com

PRODUCT SHOWCASE

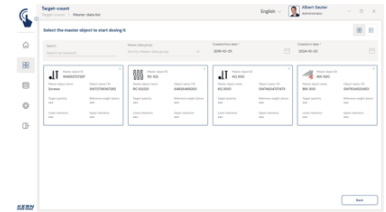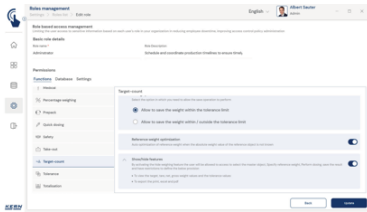

ENGLISH

DESCRIPTION

- Note: Prerequisite for this set is the basic program SET-01 Base
- SET-40 Protect is compatible with the following sets: SET-01 Base, SET-16 Batch and statistics, SET-21 Formulation, SET-32 Count, SET-33 Target count, SET-35 Assembly count and SET-31 Tolerance
- Protection of sensitive recipes: This function can be used to hide all critical recipe data in order to prevent unauthorised access to this data
- Protection of sensitive weighing data: When dosing or weighing, target weights will be shown as neutralised values by “***”, so that the result is only shown as a graphic on the basis of the target value and the lower and upper tolerance values
- Protection of sensitive object names: The name of a substance or an object will be masked for certain user groups, so that only the ID number of a substance is visible as a minimum, without the target weights
- Flexible data protection options: Depending on the user role, you can select whether only target weights should be displayed or whether the names of individual recipe components should also be displayed, for example
- Entering the total target quantity: On the startup screen, users enter the total target quantity without seeing individual target quantities or component weights
- Masked results: Target weights and actual weights are also shown as “***” in the results screen so they are masked
- Restricted display of quantity in the memory: In the dynamic memory there will be no quantity information displayed for certain user groups
- Transparent display for Admins: Clear text display of the target quantities or weights as well as the actual weights which have been weighed is possible at any time for authorised user groups such as Admins

Options

- The central data memory function Save Server (SET-10) for additional storage of all measurement data in a central, local server directory. By doing this the measurement data of all connected EasyTouch weighing

systems as well as from all installed EasyTouch functions will be stored. A particular benefit of doing this for those users with several weighing systems is that all weighing data is consolidated in just one database and you can search for individual measurement data from several balances in just one table. The Save Server data memory is also tamper-proof and cannot be changed

- **Save-Cloud:** Has the same central data memory function for all connected ET systems as Save-Data Server. The difference is that in Save-Data Cloud, the storage location is a KERN server which is accessed over the internet, instead of a server on the local network. The advantage with Save-Data Cloud is that the setup is completely automatic and you do not need a network administrator, KERN SET-101
- **EasyTouch Quick Key (SET-311):** This function enables you to reserve specific master data objects for specific user roles. The selected objects will be displayed to these users as quick-access keys on the screen

License information

- **Licensing:** One license can be operated on up to four terminal devices (PC, laptop, tablet) at the same time, working independently
- **User:** You can store as many users as you need in one license
- **Balances:** You can store and operate as many balances as you need in one license
- **Communication between balance/terminal device:** The balance(s) can communicate with the PC, laptop or tablet by serial connection, USB, Bluetooth, Ethernet or WiFi

SPECIFICATIONS

Article codes, product groups & Status information

Model article-code	SET-40
Model-Serie	ETS
Assortment	KERN
Product Type	Software
Product group	EasyTouch
Product usage type	Main item;Accessories/Services
EAN code (Model)	404576 1476547
Customs tariff number	85234910

Delivery time	1 d
---------------	-----

FUNCTIONS

Standard



Option



Technical data - Packing & Shipping