


Software EasyTouch

SET-18 Container Master

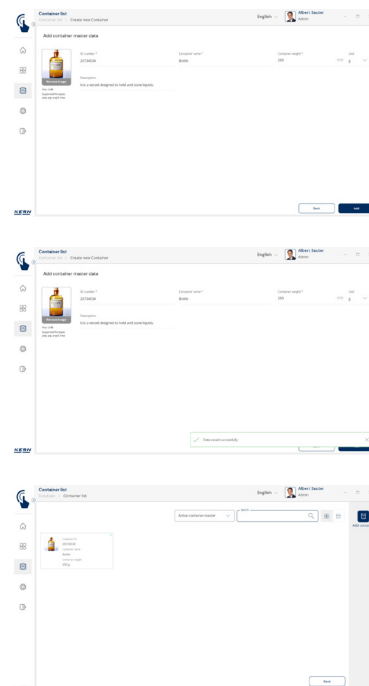

EasyTouch Container Master – Centralised container management

Features

- Note: Prerequisite for this set is the basic program SET-01 Base
- Container management: Different items are stored and weighed in similar containers. The tare weights of these containers can be managed centrally with SET-18 Container Master. If the weight of a container changes, the new tare weight is simply updated in Container Master and automatically applied to all associated items and functions in the application
- Fixed assignment: A container can be permanently assigned to each master data object to ensure consistency
- Container display: When using a master data object, in EasyTouch users see the permanently assigned container or a manually selected container. This reduces the risk of confusion and improves accuracy in the work process
- Container selection: Users can also select a container during the weighing process without having to define it beforehand, for example by scanning the container ID with the barcode. This helps to avoid incorrect tare entries

STANDARD



OPTION

