
KERN Software BalanceConnection SCD-4.0

Flexible recording or transfer of measurements, particularly to Microsoft[®] Excel

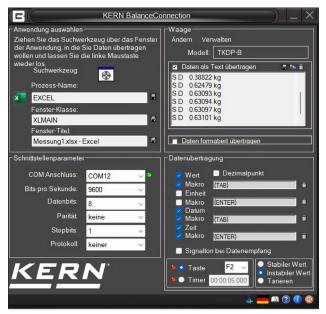

Programme start and quick configuration of the RS232/USB interface

🖲 Abschalten 📈 😝 🛻 隆 👶		
Geräte und Protokolle 🕴 🔻	Ausgebenetheden X	
🗣 Gerät hinzufügen 😿 🙀 Gerätemodelle verwalten	🛶 Ausgabe hinzufügen - 🕂 Filter - 💋 🕕 🌜	
−6.41 −6.40 −6.4	Augustantini Outperformation Outp	
	AntherMateur Greek. X	
	Andrea Maran Tana M. X Handagan (1970) Andrea	

Full overview of devices, outputs, queries, interfaces, etc. directly in the output window

CAL-INI			
TYPE	PCJ 600-2M		
SN	WD230058116		
BALID	KERN 1844		
ALIID	41		
DATE	2024 Jul 10		
TIME	14:32:19		
REF =	300.00 g		
BFR =	299.96 g		
AFT =	300.00 g		
-COMPLETE			
-SIGNATURE-			
MEMPRT			
A CONTRACTOR OF A CONTRACTOR OFTA CONTRACTOR O			
MEMPRT B MemID	42		
MEMPRT B Date	2024-07-10		
	14:32:39		
MEMPRT B Gross	233.80 g		
MEMPRT B Net	233.80 g		
MEMPRT B Tare	0.00 g		
MEMPRI B Time MEMPRI B Gross MEMPRI B Net MEMPRI B Tare MEMPRI B Range	1		
MEMPRT A			
MEMPRT			
MEMPRT B MemID	43		
MEMPRT B Date	2024-07-10		
MEMPRT B Time	14:33:01		
MEMPRT B Time MEMPRT B Gross	233.76 g		
MEMPRT B Net	222.72 g		
MEMPRT B Tare	11.04 g		
MEMPRT B Range	1		
MEMPRT A			
Send	~ I		
Text	\sim		
85 - 🗎 🧹 🔕 🔺	🔹 🖓 🗑 🎻 🔇		

Query measurement results, calibration protocol or alibi memory data via remote control commands. Can be saved as a text file.

Do you have any questions?

The team from our Industry 4.0 centre of expertise will be happy to help: e-mail: etsales@kern-sohn.com

Tip: Or download a free demo version of the BalanceConnection software: http://balanceconnection.kern-sohn.com/demo

Features

- For operating systems Windows 10, 11, -.NET Framwork 4.8.1 required
- Supports balances measuring devices with RS 232, RS 485, Bluetooth, LAN or WiFi network (TCP/UDP/IP)
- Highly flexible formatting of the output (any order, formatting and rounding), particularly recording of date and time for every value transferred, if required
- Any number of devices/interfaces can be connected, as well as simultaneous and synchronised recording of several balances
- Key-activated or time-controlled interrogation of measurements or trigger of device functions, also for continuous recording
- The interface protocols for KERN balances are already predefined (standard configuration)

- Compatible interface cable included when you order a KERN balance at the same time
- Many different transfer and recording options:
 Microsoft[®] Excel/Access/Word in the
 - foreground or background - Other Windows applications (through key simulation), e.g. shipping software or
 - ERP system (SAP, Sage, etc.)
 - File recording (e.g. as text or CSV file)
 - Print out on text or label printer
- Screen output such as large display, line chart (drying curve for moisture analyser), histogram, etc.

Benefits

- Typing errors during e.g. manual transfer are avoided
- Automatic recording, e.g. for long-term tests
- Supports GLP-compliant processes
- Can be used across different branches and applications
- Also available immediately as a download license
- Scope of supplies: Download link for 1 license, KERN SCD-4.0-DL